


NEW FPLOC INSTRUCTIONS

- 1. Place templates on force plates.
 - a. Plate 1 has the plate one square template
 - b. Plate 2 has the plate two square template
 - c. Plate 3 has the triangular frame
 - i. The corner of the triangle goes over corner 1 of the force plate
 - ii. Marker two is the marker in the clockwise direction from Marker 1
 - iii. Marker three is the other corner
 - iv. Marker four is the the marker on the cross bar of the template
- 2. Collect data using the Trial Type of FPLOC3
- 3. Close the Current Session
- 4. Minimize the Vicon 370 software
- 5. Go to NIHEXE/Cmotion/Fploc.exe
- 6. Run FPLOC.exe
- 7. Go to Properties Specify the Target Template Location

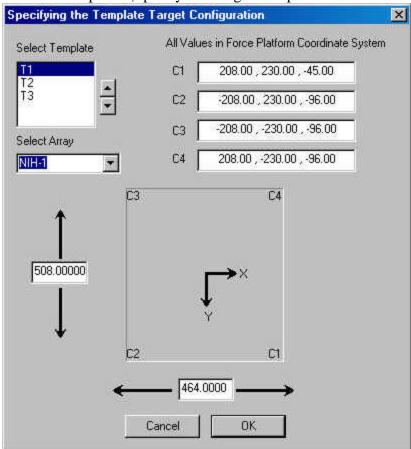


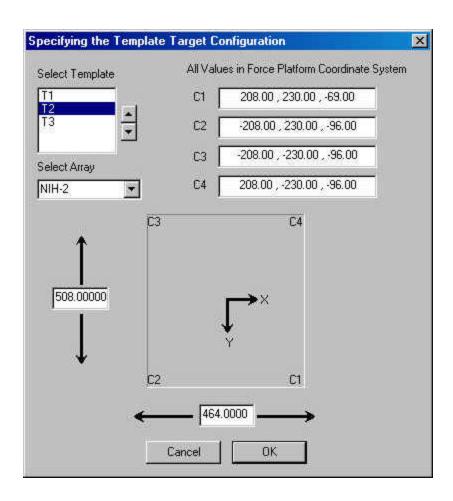

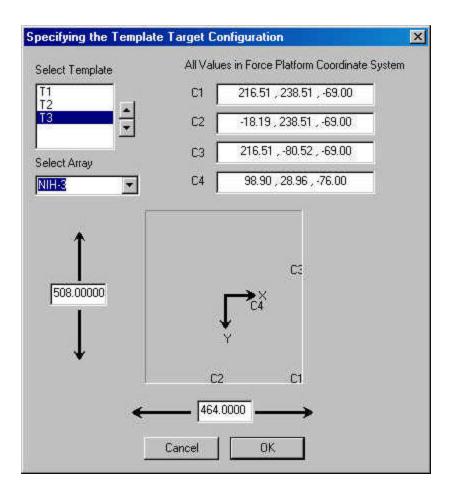


Make sure corner locations are correct. (Use Values above to check for a standard configuration with three forceplates)


- 8. Select OK
- 9. Press the open button

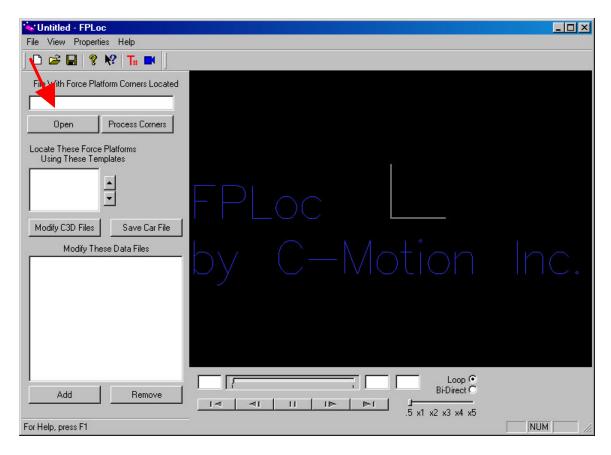

10. Select C3D file that contains the templates

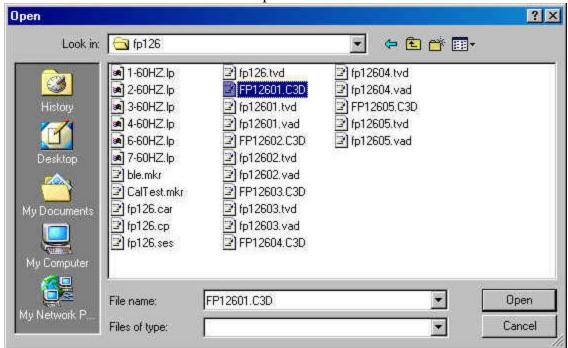


- 11. Select Open
- 12. View Template Markers in Visualization Window
- 13. Press Process Corners

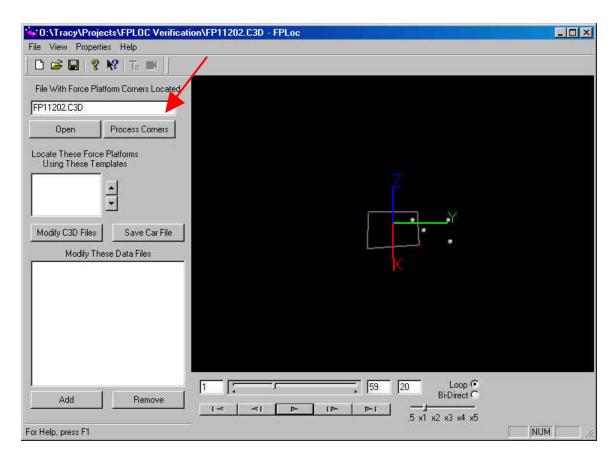

14. View Force Plates and Template Locations in Visualization Window

15. Press Save Car File

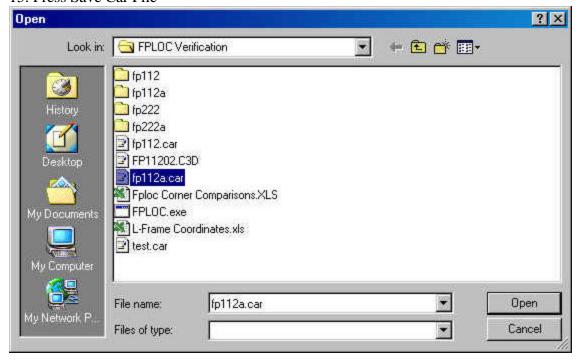

- 16. Open Car file
- 17. A window will appear that says data was saved to car file 18. Check using Vicon System that force plates were located properly